


Vigilant Centerium

Data management platform

Part of the **Agilia® Connect** Infusion System

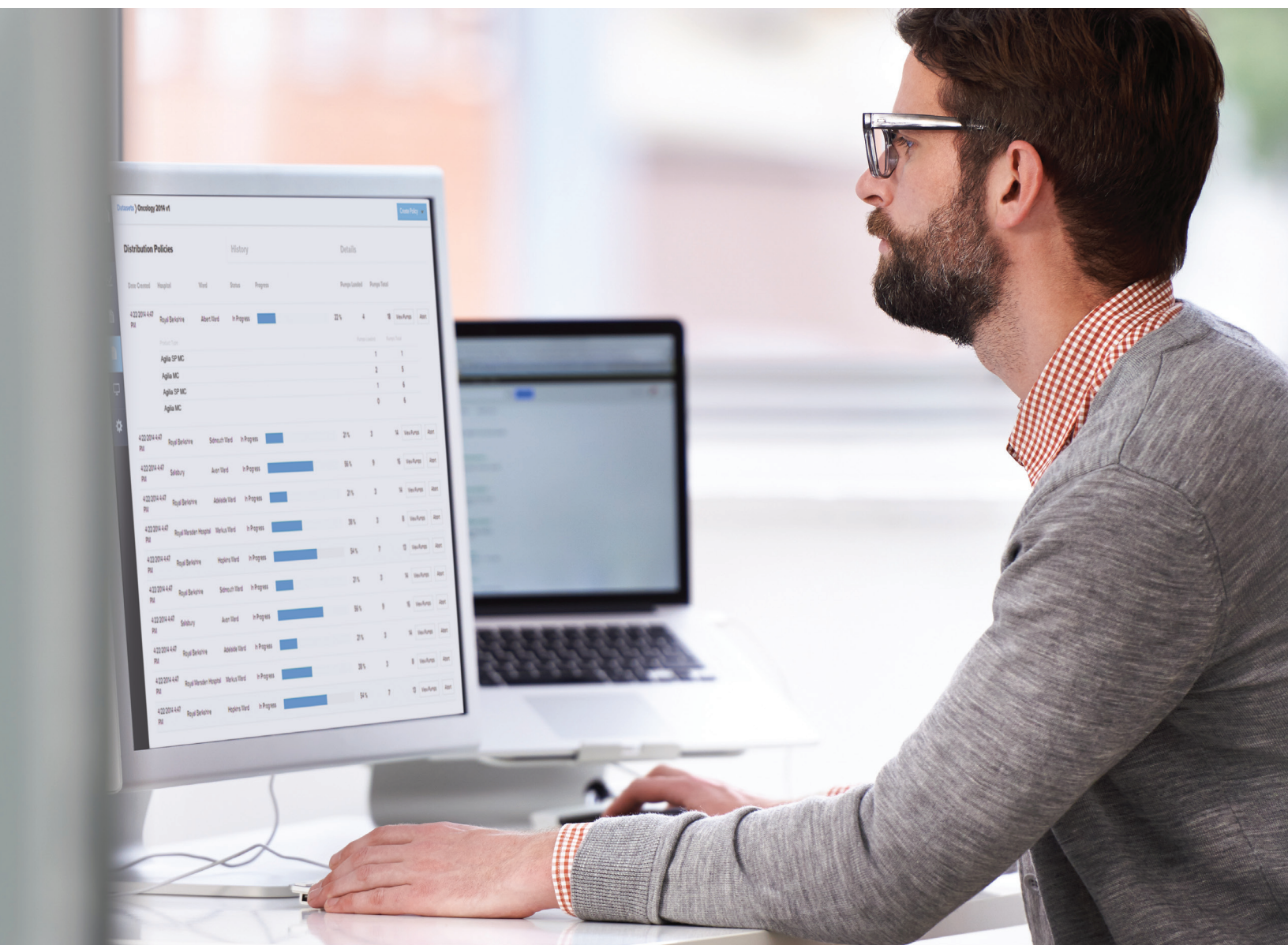

Your data management ally

The web-based Vigilant Centerium platform provides your team with wireless connectivity between the family of Agilia® Connect pumps and Vigilant Master Med and Vigilant Insight software, enabling centralized distribution of drug libraries and warehousing of infusion data for reporting and analysis.

Technically Packed with powerful cybersecurity

Built to state-of-the-art cybersecurity practices, the confidentiality and integrity of all infusion-related data that moves through Vigilant Centerium has been confirmed through an independent vulnerability test by a 3rd party testing lab.¹

Functionally Refined for intuitive operation

Vigilant Centerium's anytime, anywhere access and user-friendly interface helps ease your team's initial training and everyday use.

Practically Engineered for infrastructure efficiency

Vigilant Centerium has been designed to allow your team to maintain an optimal network load – even for 4500 pump – by fine-tuning communication frequency to best balance data availability needs and wireless bandwidth.

Exceptionally Supported to drive best practices

From initial implementation to ongoing maintenance and repair, our technical and clinical service team shows your users how to put Vigilant Centerium to best use, so they can get started quickly and work confidently.

¹. Data on file.



Dedicated Team of Experts

Account Manager
Project Manager
Clinical Pharmacy Services Liaison
Clinical Nurse Liaison
IT Service Liaison

SCP certification, 10 years running

Our technical team's commitment to annual certification by the global standard in field service quality – Service Capability & Performance (SCP) – whose audits cover a broad scope of practices, ensures installation and training, as well as ongoing service, are provided by individuals with deep knowledge.

Export all pump and infusion data in

XLS
for customized reporting

Large single-instance capacity

4,500
pumps²

Optimize infusion management

Centralized Information comes from Agilia® Connect pumps for analysis and reporting.

At-a-glance data visualizations helps to focus on pumps that may not have received the latest drug library update or communicated their most recent infusion information.

Virtualization can be enabled through VMware ESX and VMware Sphere environments.

Simple to integrate with your organization's IT systems.

Enhance fleet monitoring

Status monitoring of all installed Agilia Connect pumps.

Test, training & production environments can be created within a single instance of the software.

Implement advanced cybersecurity

Mutual Authentication using Pre-Shared Key (PSK) with TLS1.2 using AES128 (encryption) and SHA-256 (data integrity)

Controlled user access integrated with MS Active Directory.

Digital certificates for secured communications with other applications.



SCP is a registered trademark of Service Strategies.

Globally **Recognized as a proven partner**

Fresenius Kabi's medical device division is a leader throughout the world, with more than 35 years of infusion technology experience and an installed fleet of more than a million infusion pumps. Fresenius Kabi is also a steadily-growing force in the North American generic injectables space, having worked shoulder-to-shoulder with acute- and alternate-care institutions over the last decade to help ensure the availability and affordability of essential IV drugs for anesthesia and analgesia, critical care, and oncology.



**FRESENIUS
KABI**

caring for life

Fresenius Kabi, USA, LLC
Three Corporate Drive
Lake Zurich, Illinois 60047
www.fresenius-kabi.us/agiliacconnect

©2022 Fresenius Kabi USA, LLC. All Rights Reserved.
2601-ACIS-02-03/22